





## LoopMaster<sup>®</sup> - Sochman Snare

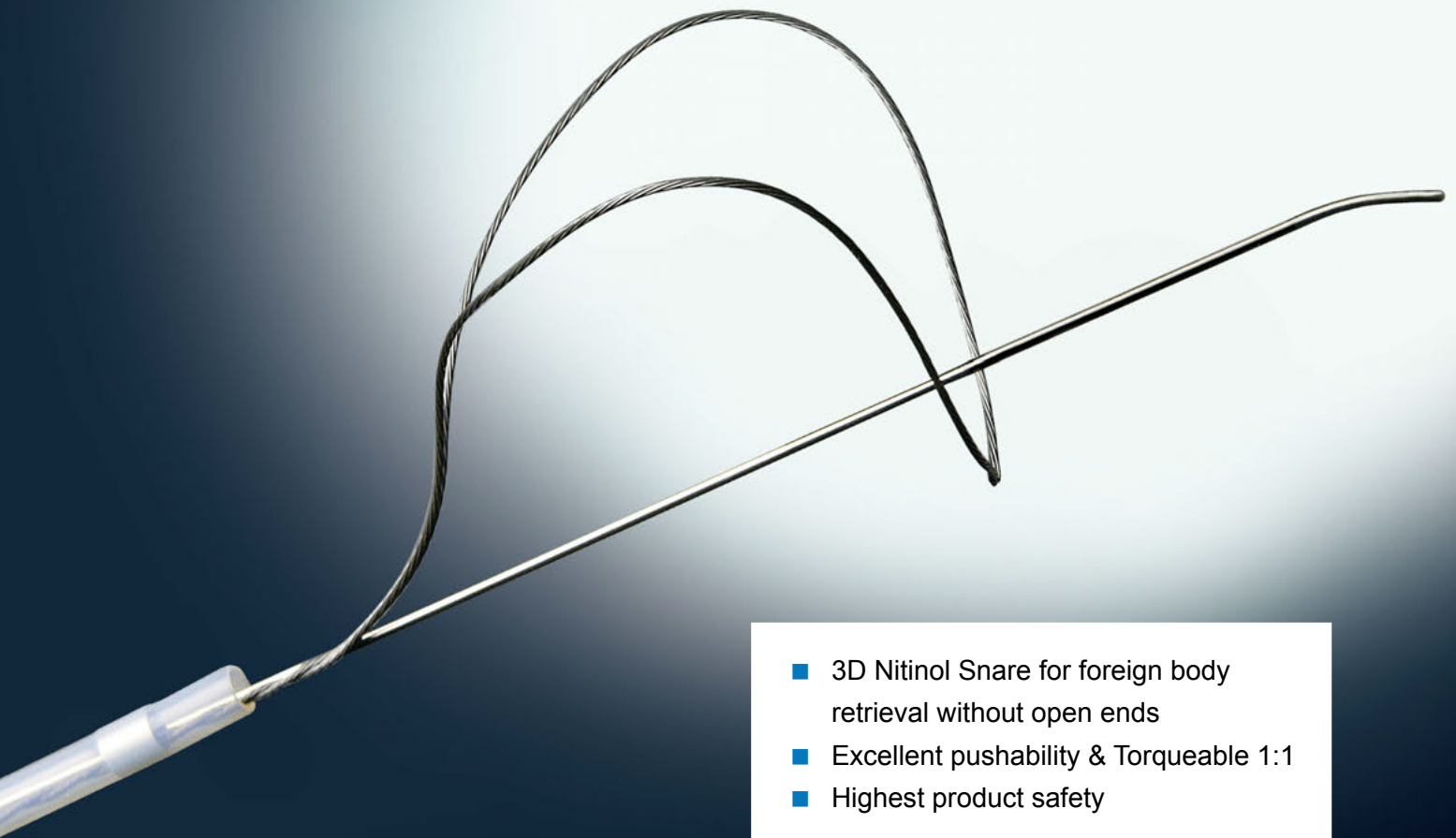

- 3D Nitinol Snare for foreign body retrieval without open ends
- Excellent pushability & Torqueable 1:1
- Highest product safety

# LoopMaster<sup>®</sup> - Sochman Snare

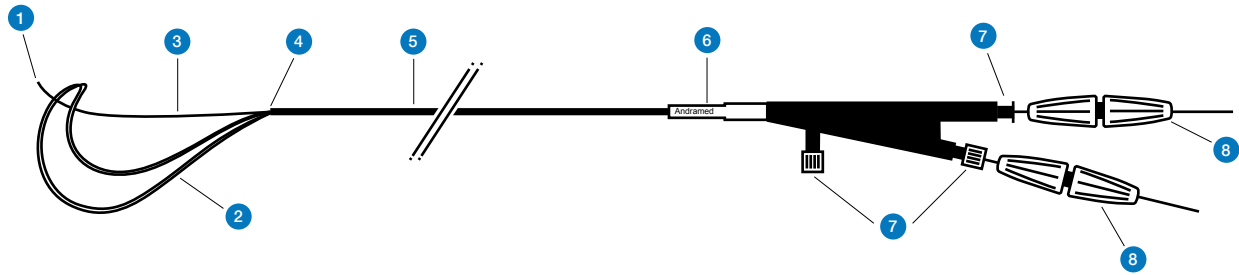


- 1 Atraumatic distal end at the catching wire
- 2 Nitinol Loop
- 3 Catching wire
- 4 X-Ray marker
- 5 Flexible shaft
- 6 Bend protection at the Luer
- 7 Luer connector
- 8 Large torquers for simple control of the catching Snare and the catching wire

## Retrieval Technique for foreign bodies



- 1 Open the Snare and place it around the foreign body to be removed.

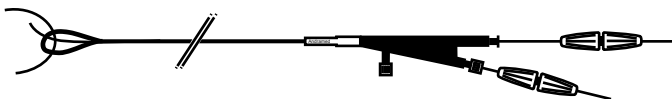

- 2 By closing the Snare, the foreign body is enclosed between the catching wire and the catching Snare.



- 3 This makes it possible to securely remove the foreign body from the body.

## Order Information LoopMaster<sup>®</sup>

Product Code	Loop- Diameter	Length (cm) Loop/Guidewire	Introducer Sheath (cm)	Introducer Sheath (F)
ASL-25	25 mm	130/150	110	8

