


Bone Biopsy Systems




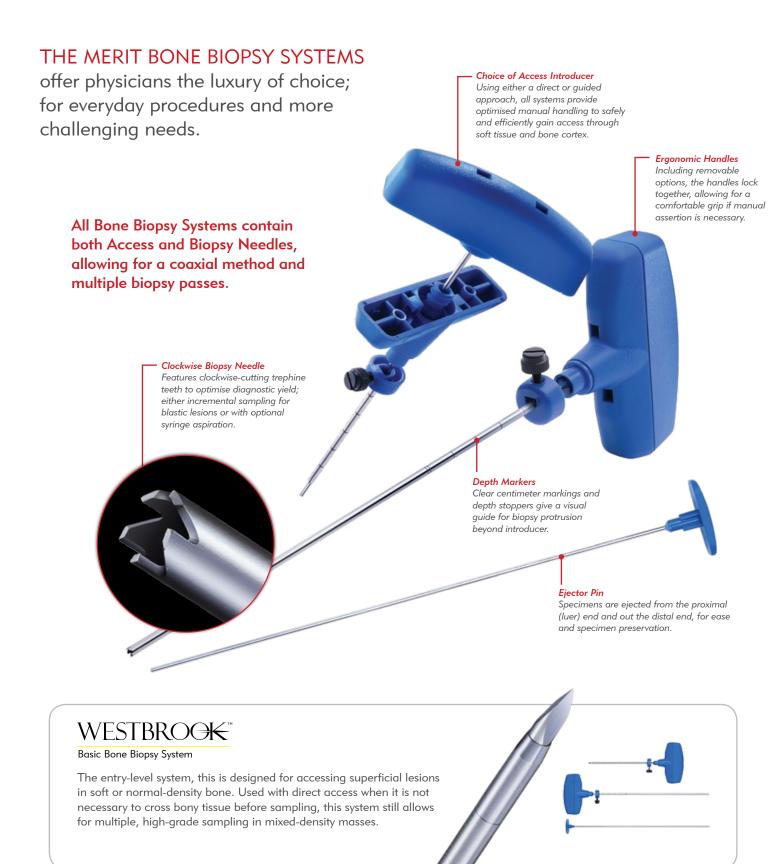

MADIS⊖N™

Perforating Bone Biopsy System


This direct access system has the patented hard bone perforating cannula and drilling technology if required to traverse hard bone and reach deep-set or small cortical abnormalities. A range of sizes make this a popular choice for both axial and appendicular cases, and the drill is locked at a set-distance whilst drilling to maximise accuracy in cannula placement and sampling precision.

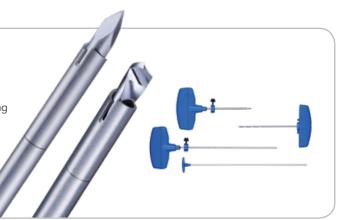

MADISON* MINI

Small Bone Biopsy System


The Madison's 'little sister', this direct-access 13G system has no drill but features a perforating cannula with removable handle. This minimises exterior weight and size and is extremely useful for small rounded bones with little soft-tissue surroundings.

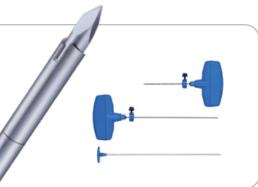

KENSINGTON*

Guide Wire Bone Biopsy System


Using guided access originally inspired by the Seldinger technique, this is the original coaxial system. With a smooth, tapered access needle placed directly over a 20G wire through the line of anaesthesia, this choice assists in creating a less traumatic pathway through soft tissues, and is designed to lesson exposure by minimising repeat scanning.

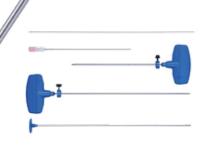

HUNTINGT⊖N™

Guided, Perforating Bone Biopsy System

Engineered to improve initial purchase on the cortex even in oblique approaches, the tri-angled drill-tip introducer provides a collection channel for bone evacuation upon first contact with bone. Accommodating the thin guide wire when Seldinger-inspired technique is desired, moving introducer counter-clockwise through soft tissue will minimise trauma, yet turning clockwise will immediately engage in cortex. An optional drill is included to work with the perforating cannula to further traverse bone where necessary.

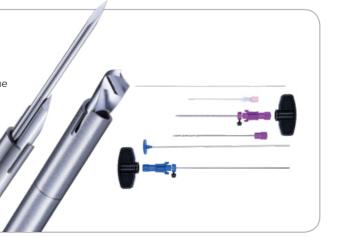

Ordering Information



Systems in 11G contain 13G* Biopsy Needles.

Systems in 13G contain 15G* Biopsy Needles.

Catalog Number	Merit Bone Biopsy Systems	Product Size*
KBD1110	Westbrook [™] Bone Biopsy System	11G, 10cm
KBD1115	Westbrook [™] Bone Biopsy System	11G, 15cm
KBC1110	Kensington [™] Guide Wire Bone Biopsy System	11G, 10cm
KBC1115	Kensington™ Guide Wire Bone Biopsy System	11G, 15cm
KDPD1165	Madison™ Perforating Bone Biopsy System with Drill	11G, 6.5cm
KDPD1185	Madison™ Perforating Bone Biopsy System with Drill	11G, 8.5cm
KDPD1110	Madison™ Perforating Bone Biopsy System with Drill	11G, 10cm
KDPD1115	Madison™ Perforating Bone Biopsy System with Drill	11G, 15cm
KDPD139	Madison™ Perforating Bone Biopsy System with Drill	13G, 9cm
KDP136	Madison™ Mini Bone Biopsy System	13G, 6cm
OM1165	Huntington [™] Guided, Perforating Bone Biopsy System with Drill	11G, 6.5cm
ОМ1110	Huntington [™] Guided, Perforating Bone Biopsy System with Drill	11G, 10cm
ОМ1115	Huntington™ Guided, Perforating Bone Biopsy System with Drill	11G, 15cm

Catalog Number	Preston [™] Bone Accessories	Product Size*
ANP1165	Preston™ Hard Bone Access Needle with Drill	11G, 6.5cm
ANP1185	Preston™ Hard Bone Access Needle with Drill	11G, 8.5cm
ANP1110	Preston™ Hard Bone Access Needle with Drill	11G, 10cm
ANP1115	Preston™ Hard Bone Access Needle with Drill	11G, 15cm
KTP1110	Preston [™] Hard Bone Introducer	11G, 10cm
KTP1115	Preston [™] Hard Bone Introducer	11G, 15cm
CD1110	Preston™ Bone Biopsy Needle	13G, 20cm
CD1115	Preston [™] Bone Biopsy Needle	13G, 24cm
ML1110	Preston [™] Long Drill	13G, 18.6cm
ML1115	Preston [™] Long Drill	13G, 23.6cm
TV1110	Preston™ Transpedicular Introducer	11G, 10cm
TV813	Preston™ Transpedicular Introducer	8G, 13cm
MV1110	Preston™ Transpedicular Drill	13G, 13.7cm
CR1320	Preston [™] Bone Fillers	13G, 20cm
CR1020	Preston [™] Bone Fillers	10G, 20cm

*Product size refers to access needles when describing Merit Bone Biopsy Systems. Gauge sizes may be closest approximate.

Before using refer to Instructions for Use for indications, contraindications, warnings, precautions, and directions for use.



Understand. Innovate. Deliver.™

Merit Medical Systems, Inc. 1600 West Merit Parkway South Jordan, Utah 84095 1.801.253.1600 1.800.35.MERIT

merit.com

Merit Medical Europe, Middle East, & Africa (EMEA) Amerikalaan 42, 6199 AE Maastricht-Airport

The Netherlands +31 43 358 82 22

Merit Medical Ireland Ltd.

Parkmore Business Park West Galway, Ireland +353 (0) 91 703 733

Australia

+61 (0) 3 8373 4943

Austria

0800 295 374

0800 72 906 (Dutch) 0800 73 172 (Français)

+55 11 4561 - 0062

Canada 1.800.364.4370

China +86-10 8561 0788

Denmark 80 88 00 24

Finland 0800 770 586

0800 91 60 30

Germany 0800 182 0871

Hong Kong +852-2207 6888

India

+91-80-41223376 Ireland (Republic) 1800 553 163

Italy 800 897 005 8002 25 22

Netherlands 0800 022 81 84

Norway 800 11629

Portugal 308 801 034

Russia

+7 495 221 89 02 Singapore +65 67536311

+8210 2205 0707

Spain +34 911238406

Sweden

020 792 445 Switzerland

(Deutsch) +41 22 518 02 30

(Français) +41 22 518 02 52

(Italiano) +41 22 518 00 35

UK 0800 973 115

403991101/001 ID 010218